

Actitouch Delivers Accurately Targeted and More Consistent Compression versus Multi-Layer Bandages

Original Article: Comparison of 4-layer bandages and an Adaptive Compression Therapy device on intended pressure delivery

Journal: *Journal of Wound, Ostomy, & Continence Nursing*

Publication Date: September/October 2015

Authors: Harvey N. Mayrovitz*, Hugo Partsch†, Wolfgang Vanscheidt‡

INTRODUCTION

Sustained compression provided by compression bandaging is the underpinning for treatment of chronic venous insufficiency and venous leg ulcers (VLU)¹⁻⁵, but is challenged by a loss of pressure over time due to a decrease in limb volume in both healthy and edematous individuals⁶. The commonly used “four-layer bandage” (4LB) is designed to provide a pressure profile of 40 mmHg at the ankle graduating to 17 mmHg at the upper calf when applied correctly per manufacturer’s instructions. Actitouch® Adaptive Compression Therapy is a 4-chamber pneumatic compression device which offers both intermittent and sustained compression⁷. Pressures are monitored by a microprocessor which measures and adjusts inflation levels every 30 minutes to achieve targeted interface pressures of 40, 30 and 20 mmHg distal to proximal.

OBJECTIVE

To compare interface pressure profiles of Actitouch therapy with traditional 4LB and compare the ability of each to sustain graduated compression over time.


METHODS

Actitouch and 4LBs were evaluated for interface pressure magnitude, variability, and consistency over time on 12 healthy individuals (24 legs). Subjects wore Actitouch on one leg while a new 4LB was applied daily on the opposite foot by a single, experienced healthcare professional. A PicoPress® Compression Measurement System measured interface pressures at three sites on the lower leg (Figure 1). Pressures were recorded with subjects seated and standing at time of device application and at 1, 4 and 8 hours of device wear.

RESULTS (Figure 2)

At 1 hour, the 4LB showed a significant ($P < 0.001$) decline in pressures while seated at all three pressure sites while Actitouch pressures showed no decline. By hour 8 of device

Figure 1

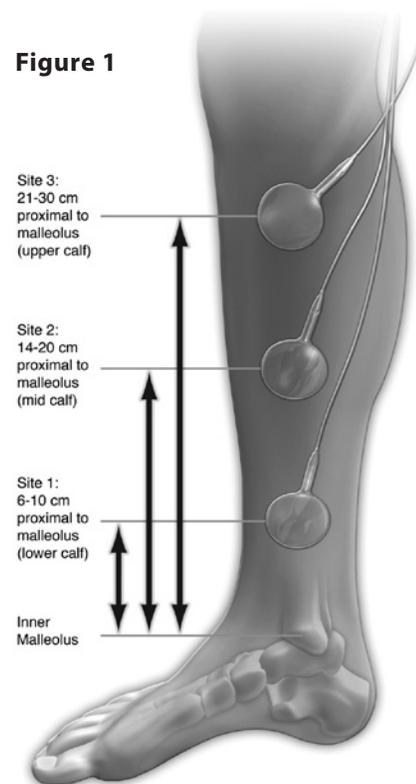


wear, the largest average seated pressure decrease for the 4LB occurred at site 2 and was 9.1 mmHg, while Actitouch showed an increase of 4.9 mmHg at the same site.

DISCUSSION

Static and dynamic (also called “working”) pressures are beneficial in the treatment of VLUs⁸⁻¹¹. Static compression is most often applied using compression bandaging, with dynamic pressure produced by the calf muscle working against the counterforce of this bandaging when standing or walking. Intermittent pneumatic compression (IPC) simulates this working pressure in patients with low activity levels or an insufficient calf muscle pump¹²⁻¹⁵. While both static and dynamic compression are important, most devices

Figure 2 – Seated Pressures

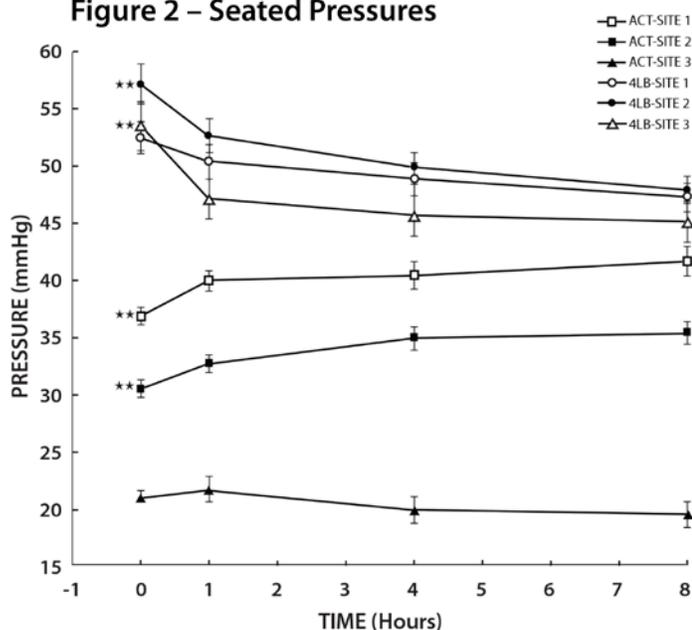

are developed to apply only one type of pressure. This data shows Actitouch is able to provide targeted static compression levels in addition to dynamic pressure using IPC with one device regardless of patient ambulation.

CONCLUSION

Actitouch delivers a more consistent and predictable sustained interface pressure when compared to a 4LB. The 4LB demonstrated significant pressure loss.

KEY POINTS

- Even with daily application, a 4LB loses interface pressure immediately and continuously – a mean of 16% at 8 hours
- Actitouch delivers consistent target interface pressures immediately and over time
- Both static and dynamic compression are important in treatment of VLU
- Actitouch may be self-applied by a patient and achieves consistent pressures

About the Authors

*Harvey N. Mayrovitz, PhD, Professor of Physiology, College of Medical Sciences, Nova Southeastern University, Ft Lauderdale, FL
 †Hugo Partsch, MD Emeritus Professor of Dermatology, Department of Dermatology, Medical University of Vienna, Vienna Austria
 ‡Wolfgang Vanscheidt, MD Professor of Dermatology, Facharzt für Dermatologie Phlebologie Allergologie, Freiburg, Germany

References

- The full article can be found in the *J Wound Ostomy Continence Nurs.* 2015;42(5):1-6 Published by Lippincott Williams & Wilkins.
1. O'Meara S, Cullum N, Nelson EA, Dumville JC. Compression for venous leg ulcers. *Cochrane database of systematic reviews.* 2012;11:CD000265.Pub3.
 2. Nikolovska S, Arsovski A, Damevska K, Gocev G, Pavlova L. Evaluation of two different intermittent pneumatic compression cycle settings in the healing of venous ulcers: a randomized trial. *Med Sci Monit.* Jul 2005;11(7):CR337-343.
 3. Moffatt CJ, Edwards L, Collier M, et al. A randomised controlled 8-week crossover clinical evaluation of the 3M Coban 2 Layer Compression System versus Profore to evaluate the product performance in patients with venous leg ulcers. *Int Wound J.* Jun 2008;5(2):267-279.
 4. O'Donnell, TF, Passman MA, Marston WA, et al. Management of venous leg ulcers: clinical practice guidelines of the Society for Vascular Surgery (R) and the American Venous Forum. *J Vasc Surg.* 2014. 60(2 Suppl): p. 3s-59s.
 5. Hafner J, Botonakis I, Burg G. A comparison of multilayer bandage systems during rest, exercise, and over 2 days of wear time. *Arch Dermatol.* Jul 2000;136(7):857-863.
 6. Damstra RJ, Brouwer ER, Partsch H. Controlled, comparative study of relation between volume changes and interface pressure under short-stretch bandages in leg lymphedema patients. *Dermatol Surg.* Jun 2008;34(6):773-8.
 7. Harding KG, Vanscheidt W, Partsch H, Caprini JA, Comerota AJ. Adaptive compression therapy for venous leg ulcers: a clinically effective, patient-centred approach. *Int Wound J.* May 7 2014.doi: 10.1111/iwj.12292.
 8. Partsch H. Venous narrowing by compression of the lower extremities: a prerequisite for improving venous hemodynamics. *Vasa.* 2014 Jul;43(4):235-7.
 9. Mayrovitz HN, Delgado M, Smith J. Compression bandaging effects on lower extremity peripheral and sub-bandage skin blood perfusion. *Ostomy Wound Manage.* Mar 1998;44(3):56-60, 64 passim.
 10. Mayrovitz HN, Sims N. Effects of ankle-to-knee external pressures on skin blood perfusion under and distal to compression. *Adv Skin Wound Care.* Jul-Aug 2003;16(4):198-202.
 11. Partsch H, Clark M, Mosti G, et al. Classification of compression bandages: practical aspects. *Dermatol Surg.* May 2008.
 12. Warwick D, Dewbury K, Forrester A. Intermittent pneumatic compression. A comparison of femoral vein velocity with five different devices. *Int Angiol.* Aug 2013;32(4):404-409.
 13. Iwama H, Obara S, Ohmizo H. Changes in femoral vein blood flow velocity by intermittent pneumatic compression: calf compression device versus plantar-calf sequential compression device. *J Anesth.* 2004;18(3):232-233.
 14. Delis KT, Slimani G, Hafez HM, Nicolaidis AN. Enhancing venous outflow in the lower limb with intermittent pneumatic compression. A comparative haemodynamic analysis on the effect of foot vs. calf vs. foot and calf compression. *Eur J Vasc Endovasc Surg.* Mar 2000;19(3):250-260.
 15. Janssen H, Trevino C, Williams D. Hemodynamic alterations in venous blood flow produced by external pneumatic compression. *J Cardiovasc Surg (Torino).* Oct 1993;34(5):441-447.
 16. Partsch H, Mosti G. Pressure-time integral of elastic versus inelastic bandages: Practical implications. *EWMA Journal* 2013;13(2):15-17.

Tactile Medical

1331 Tyler Street NE, Suite 200
 Minneapolis, MN 55413 USA
 T: 612.355.5100
 F: 612.355.5101

Customer Service

Toll-Free Tel: 866.435.3948
 Toll-Free Fax: 866.435.3949
 Hours: Monday through Friday, 7 a.m. – 7 p.m. CT
www.tactilemedical.com

Tactile
MEDICAL™
 HEALING RIGHT AT HOME